BORT Air Walker kurz

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Med
Mitten im
gesunden
Leben

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

REF 100 350

Illustrationen
Figures


1

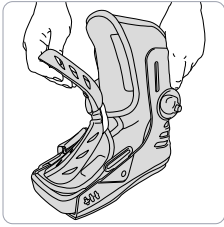
2

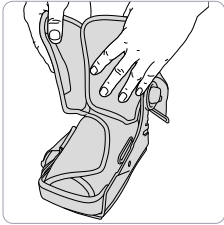
3

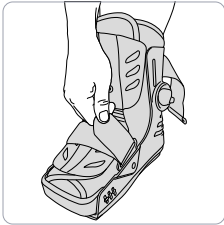4

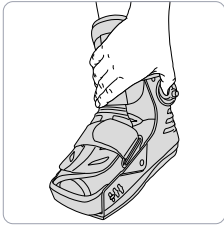5

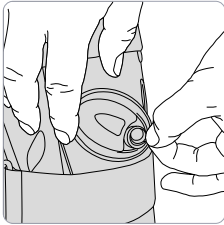Sprachen Languages

DE	deutsch	Gebrauchsanweisung	04
EN	english	Instructions for use	07
FR	français	Mode d'emploi	10
ES	español	Instrucciones de uso	13
IT	italiano	Istruzioni per l'uso	16
NL	nederlands	Gebruiksaanwijzing	19
CS	český	Návod k použití	22
ET	eesti	Kasutusjuhend	25
PL	polski	Instrukcja użytkowania	28
RO	româneșc	Instrucțiuni de utilizare	31

 PDF: ga.bort.com

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Many thanks for placing your trust in a medical device from BORT GmbH. Please read the existing instructions for use carefully. If you have any questions, please refer to your physician or the specialist retailer from whom you purchased this medical device.

Intended purpose

This medical device is a lower leg foot brace for complete immobilisation of the lower leg and foot region in a specified position.

Indications

Pre-operative or post-operative/post-traumatic immobilisation of the lower leg and metatarsus region, metatarsal fracture, soft tissue injuries, distortions.

Contraindications

Thrombosis risk, extreme varicosis, peripheral arterial occlusive disease (PAOD), lymph drainage disorders, also unclear soft tissue swellings distal to the aid positioned, sensory and circulatory disorders of the region of the body treated, skin diseases in the part of the body treated.

Application risks/Important notes ⚠

This medical device is a prescribable product. Discuss the use and duration with your treating physician. The expert staff from whom you have received the medical device will select the appropriate size and instruct you regarding its use.

- remove the medical device prior to radiological examinations
- the Walker is appropriate for a maximum body weight of 130 kg
- if wearing at night is prescribed, avoid negative influence on the circulatory system
- in case of numbness, loosen the medical device and remove it if necessary
- in case of persistent complaints, consult the physician or specialist retailer
- use the medical device in accordance with therapeutic needs
- only use other products simultaneously after consultation with your physician
- do not make any changes to the product
- do not wear it on open wounds
- do not use in case of intolerance of one of the materials used
- no re-use – this medical aid is intended for treating one patient
- whilst wearing the brace, please neither use any creams nor ointments on or around it as they can destroy the material

Fitting instructions

Before going through steps 1–5, pull one of the two pairs of understockings provided over your foot/lower leg.

Information:

The opaque, see-through protective caps are only used for transporting the Walker. They must be removed prior to use. When ready, the protective caps no longer need be attached to the blue pumps at the side. They could well prove to be an obstacle when walking.

- 1 Open the circumferential Velcro fasteners and remove the front cover.
- 2 Open the black padding liner and place the foot and lower leg in the Walker. Pay attention to good heel positioning. Then, the black padding liner is closed, pay attention it sitting crease-free.
- 3 Now, adapt the grey front cover by attaching it centrally to the black padding liner again using Velcro. Close the 2 circumferential Velcro fasteners. Start with the fastener over the instep area and finally the lower leg fastener. If necessary, re-tighten again.
- 4 Now, fill the integrated air chambers using the blue pumps at the side. The air chambers can be filled individually. Start with the inner (medial) pump, then the outer (lateral) pump.
- 5 If you wish to reduce the air pressure, you can release air. For this, fold the grey protective cap next to the pumps upwards and press the silver-coloured button.

Accessories available:

Ref. 100 370 Understocking

Ref. 100 380 Protective cover against wetness

Ref. 100 390 Protective cover for night-time

Removal

If you wish to remove the Walker, fold the grey protective cap next to the pumps upwards and press the silver-coloured button. Then, open the circumferential Velcro fasteners and remove the front cover. Finally, open the black padding liner and remove your foot from the Walker.

Tip: Close the Velcro fastener in order to avoid early wear to the Velcro surface.

Material composition

Polyester (PES), polyurethane (PUR), polyamide (PA)

The sewn in textile label provides the precise material composition.

Cleaning information

🧼 Hand-washing ☒ Do not clean chemically ☒ Do not bleach ☒ Do not dry in a tumble dryer

☒ Do not iron

Hand wash the shoe using lukewarm water and mild detergents.

Guarantee

The legal regulations of the country in which you acquired the product apply to the purchased product. Please contact your specialist retailer if you suspect a warranty claim. Please clean the product before submitting a warranty claim. If the enclosed instructions for use have not been properly observed, the warranty may be impaired or cancelled. The warranty does not cover use of the product inappropriate for the indication, non-observance of application risks, instructions and unauthorised modifications to the product.

Useful life/Lifetime of the product

The lifetime of the medical device is determined by its natural wear and tear if treated appropriately and as recommended.

Duty of notification

If a serious worsening of the state of health occurs when using the medical device, you are obliged to notify the specialist dealer or us as the manufacturer and the Medicines & Healthcare products Regulatory Agency (MHRA).

You can find our contact information in these instructions for use. You can find the contact information for the appointed authority for your country under the following link: www.bort.com/md-eu-contact.

Disposal

Upon the termination of use, the product must be disposed of in accordance with the corresponding local requirements.

Declaration of conformity

We confirm that this device conforms with the requirements of REGULATION (EU) 2017/745 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL. You can find the current declaration of conformity under the following link: www.bort.com/conformity

Status: 10.2019

Medical device |  Single patient – multiple use

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